





BMI Scale

User Manual

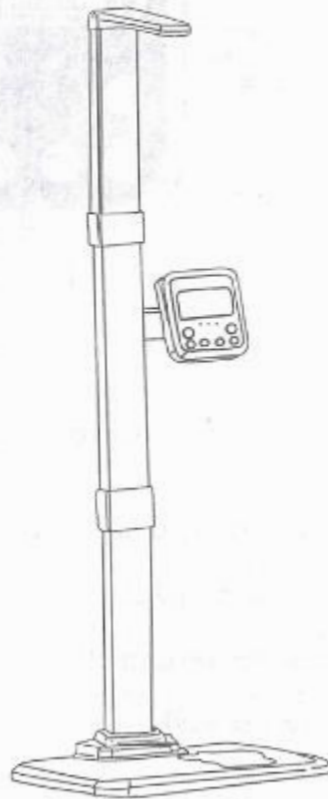

Customer Service: (888) 439-3410
houseables.com


Display Guide

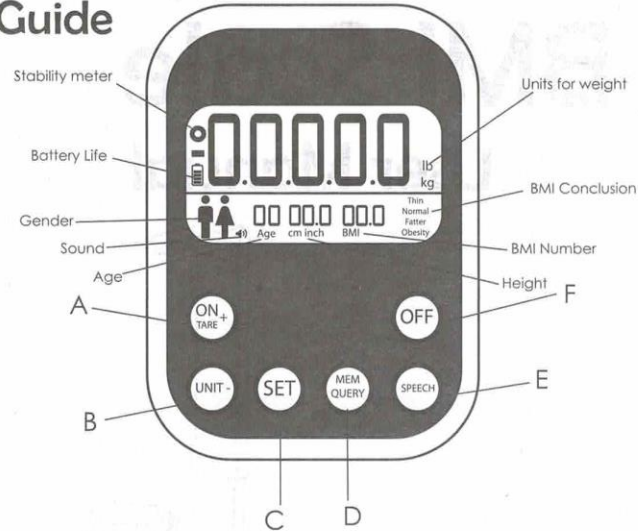

- A.** Turns the display on and increases personal parameters
- B.** Changes unit conversion and decreases personal parameters
- C.** Confirm personal parameters
- D.** Saves queries and stats
- E.** Turns the voice on and off
- F.** Turns the display off


Scale Guide

- A.** Top column with sensor
- B.** Middle column with display
- C.** Bottom column
- D.** Platform
- b1.** Rubber sleeve
- c1.** Decorative shell
- E.** M5 Phillips head large screws
- F.** M6 Hexagon screws
- G.** Flange
- H.** Display

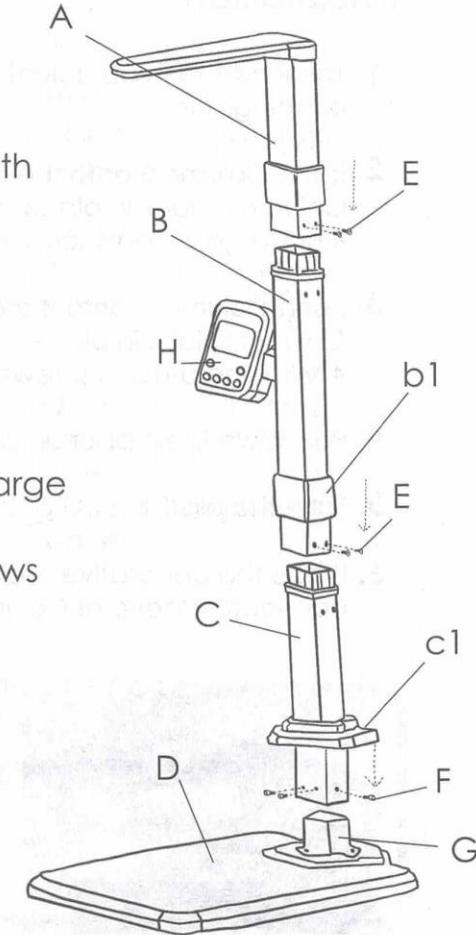

Installation

1. Take out the three columns (A, B, C), lay them on the ground
2. Insert Column B onto the connecting sleeve of Column C, lock in place by screwing in 4 M5 large b-head screws (E)
3. Insert Column A onto the connecting sleeve of Column B, lock in place by screwing in 4 M5 large b-head screws (E)
4. Pull down the rubber sleeve (b1) until it locks in place.
5. Take the platform (A) out and lay it on a flat surface
6. Place the decorative shell (c1) onto the bottom connector sleeve of Column C

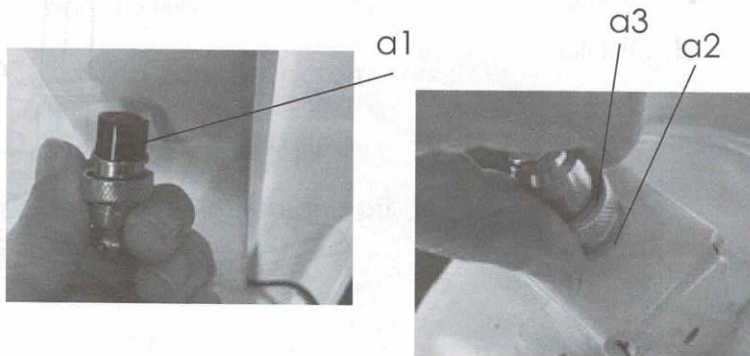

Installation

7. Connect the Aviation Plug (a1) to the Aviation Socket (a2) that is on the platform (A). Using the metal ring (a3) around the Aviation Socket (a2), screw in CLOCKWISE to lock the socket and plug together.
8. Connect the fully assembled Column onto the flange (a4) of the platform, and lock in place using 4 M6 hexagon screws (F)
9. Slide the Protective Shell (c1) downward into the corresponding assembly hole on the platform (A).
10. Lift item off of the ground to stand straight and place on a level, dry, flat space.

Note: The Display (H) can be rotated 90° by taking out the knurled screw, rotating the Display (H) 90°, and screwing the original knurled screw back into place to lock in.

Customer Service: (888) 439-3410
houseables.com

Customer Service: (888) 439-3410
houseables.com



How to Use

When using this product, ensure that you have placed it in a level, flat, and dry space. Turn it on by flipping the switch on the back of the Display (H).

Using without Personal Parameters

1. Turn the scale on using the button on the back of the display.
2. Press 'ON/Tare/+' to zero out the display
3. Stand in the center of the scale, keeping your head up, back straight, looking straight ahead. Wait for the scale to finish the height and weight measurement
4. Measurement is completed. Step off the scale and press 'OFF' to turn the display off.

Using with Personal Parameters

1. Turn the scale on using the button on the back of the display.
2. Press 'ON/Tare/+' to zero out the display
3. Press 'SET', the male/female in the bottom left of the display will begin to flash. Use 'ON/Tare/+' to select gender, press 'SET' to confirm and move on.
4. The age display should begin to flash, use 'UNIT/-' to decrease the age OR 'ON/Tare/+' to increase the age.
5. Press 'SET' to confirm and move on.

Customer Service: (888) 439-3410
houseables.com



To Save Your Data

When your measurement is finished on the scale, stay standing on it and press the 'MEM/QUERY' button to save the data

The Display should display SAVE

To Change Languages

To change the language of your scale, press and hold the 'SPEECH' button until an -E- shows. This is for English. Continue until you find your preferred language.

When the voice assistant is turned off, you will be able to hear the ultrasonic height measurer make a small static sound, when that stops it means it's done measuring

Thank you!

By purchasing this BMI scale, you have personally supported a small business in sunny Southern California! Thank you so much!

Customer Service: (888) 439-3410
houseables.com



How to Use

When using this product, ensure that you have placed it in a level, flat, and dry space. Turn it on by flipping the switch on the back of the Display (H).

Using without Personal Parameters

1. Turn the scale on using the button on the back of the display.
2. Press 'ON/Tare/+' to zero out the display
3. Stand in the center of the scale, keeping your head up, back straight, looking straight ahead. Wait for the scale to finish the height and weight measurement
4. Measurement is completed. Step off the scale and press 'OFF' to turn the display off.

Using with Personal Parameters

1. Turn the scale on using the button on the back of the display.
2. Press 'ON/Tare/+' to zero out the display
3. Press 'SET', the male/female in the bottom left of the display will begin to flash. Use 'ON/Tare/+' to select gender, press 'SET' to confirm and move on.
4. The age display should begin to flash, use 'UNIT/-' to decrease the age OR 'ON/Tare/+' to increase the age.
5. Press 'SET' to confirm and move on.

Customer Service: (888) 439-3410
houseables.com



To Save Your Data

When your measurement is finished on the scale, stay standing on it and press the 'MEM/QUERY' button to save the data

The Display should display SAVE

To Change Languages

To change the language of your scale, press and hold the 'SPEECH' button until an -E- shows. This is for English. Continue until you find your preferred language.

When the voice assistant is turned off, you will be able to hear the ultrasonic height measurer make a small static sound, when that stops it means it's done measuring

Thank you!

By purchasing this BMI scale, you have personally supported a small business in sunny Southern California! Thank you so much!

Customer Service: (888) 439-3410
houseables.com